

CRYOTHERAPY
FREEZING TECHNOLOGY
AS USED BY DOCTORS

CRYOTAG

SKIN TAG REMOVER

**FAST
EFFECTIVE
SAFE**



**READ ALL INSTRUCTIONS
BEFORE USING CRYOTAG.**

CRYOTAG is the home-use skin tag remover that allows you to remove irritating or unsightly skin tags quickly and effectively, leaving you with healthy looking skin.

FAST

Removes skin tags in as little as 1 treatment*. Most skin tags will fall off after 2 weeks.

EFFECTIVE

Clinically proven**. With its precision tip applicator and tweezers, Cryotag instantly freezes the skin tag to the core. Contains 12 applicators (12 treatments).

SAFE

Cryotag uses the same freezing technology used by doctors and dermatologists. Thanks to its precision tip applicator, Cryotag targets only the skin tag itself, leaving the surrounding skin untouched.

Cryotag is intended for the removal of skin tags ONLY. Please read Section A on "What is a Skin Tag" to help you correctly identify your growth. If you have a skin tag, read the rest of these instructions carefully. See How to Use in section H.

DO NOT USE CRYOTAG IF:

- You have any type of skin cancer or if you ever had skin cancer.
- You have or had any skin lesion diagnosed as pre-cancerous anywhere on your body.

If you treat a skin condition that is not a skin tag, you delay proper treatment. This delay may lead to a more serious skin condition or disease, and could damage healthy skin or nerves.

*Some skin tags may require additional treatments. You can treat the skin tag again, up to 4 treatment cycles.

** 1993 and 2007 clinical trials.

A. WHAT IS A SKIN TAG?

A skin tag is a small, soft, piece of skin that may be raised from the skin and hang off the skin by a connecting stalk. They are usually similar in color to the surrounding skin.

They are usually found on the neck, chest, back, under the breasts, under the arms or in the groin area. They are most common in middle aged and older adults, diabetics and overweight individuals. Skin tags are harmless and are not contagious. They usually don't cause any pain, however they can be unsightly or irritating and can cause discomfort especially if anything rubs on them, such as clothing or jewellery.

If you are unsure about whether you have a skin tag, please consult your doctor.

**Read all WARNINGS before using this product .
Only use this product if your skin tag has ALL 4 characteristics shown in the pictures below.**

Location

Skin tags MUST be located in one of the 5 locations. Only treat skin tags if they are located on the neck, chest, underarms, back or groin.



Colour

Skin tags MUST be similar in colour to the surrounding skin.



Shape & feel

Skin tags MUST be soft and flexible. You should be able to move them back and forth with your finger. They must also be raised above the skin on a stalk or stem like piece of skin.



Size

Skin tags MUST be 2mm to 6mm in size.



If your skin tag does not have ALL 4 characteristics, do not use this product. See a doctor in case you have a more serious skin condition.


For questions about the treatment of skin tags not covered in these instructions, please seek the advice of a GP or dermatologist.

B. CRYOTAG KIT CONTENT:

Cryotag contains 4 elements plus the instruction leaflet.

1. A pressurized aerosol canister with a safety shield.
2. A cap that holds the applicator.
3. 12 precision-tip applicators for 12 treatments.
4. Disposable tweezers.

This instruction leaflet contains detailed information about Cryotag, when to use it, how it works and how it can be used to remove skin tags.

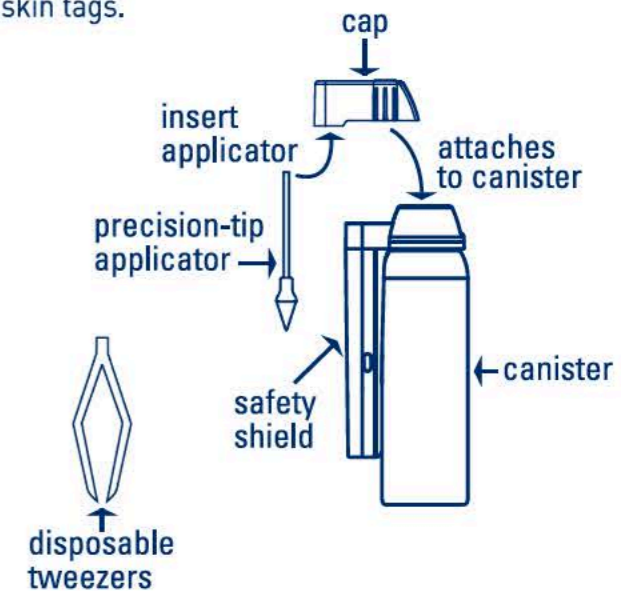

C. HOW DOES CRYOTAG WORK?

Cryotag works by instantly freezing the core of the skin tag, using the same treatment method as doctors. Freezing the skin tag destroys the tissue, which will usually cause the skin tag to fall off with 1 treatment.

- When the applicator is applied to the skin tag, the skin will turn white as a sign that it is being frozen. The treatment will feel very cold and a slight aching, stinging sensation might occur as a result of the freezing. These feelings are temporary, and should begin decreasing after the applicator is removed and as time passes. Due to the freezing, a small blister may form (sometimes containing blood) within a few days. This can cause minor discomfort but this again will disappear as time passes.

After treatment, the skin tag may darken as a sign the tissue has been destroyed.

- Most skin tags will fall off within 2 weeks of the treatment, leaving you with healthy looking skin.

D. WHEN NOT TO USE CRYOTAG

Do not use if you:

- **Are a diabetic, have poor blood circulation, have been diagnosed with blood conditions affected by extreme cold or have a weakened immune system.** Cryotag works by freezing the skin tag. People with these conditions will find that the treated area may not heal as well or completely.
- **Have a bleeding problem.** There is a risk of bleeding after the treatment.
- **Are pregnant or breast-feeding.** It is uncertain how Cryotag may affect patients under these conditions.
- **Are under the age of 18 years old.** It is uncertain how this product may affect you.
- **Do not use on skin growths that are not skin tags,** such as moles, birthmarks, dark-coloured, hairy or any unusual-looking skin conditions. If you treat a skin condition that is not a skin tag, you delay proper treatment. This delay may lead to a more serious skin condition or disease, and could damage healthy skin or nerves.
- **Do not use inside your mouth, nose, anus, genitals, buttocks or lips or on skin tags in inaccessible or hard to see places.** This may cause burns and permanent scarring of the skin.
- **Do not use on irritated or infected skin with signs of redness, itching or swelling.** It may worsen your current condition and make the infection more severe.
- **Do not use on eyes, eyelids, face or genitals.** It may cause burns, permanent scarring or blindness.

E. CAUTIONS:

- **Do NOT use an applicator other than those supplied with the product.** Always use one of the applicators supplied with the product when applying Cryotag. 1 applicator for each skin tag.

- **Do NOT treat more than 1 skin tag per location on the body at one time.** Any remaining skin tags can be treated separately with a 2 week interval between treatments. This minimizes any chance of "over freezing," which can result in burns and permanent scarring of healthy tissue.
- **Do NOT treat skin tag for longer than 40 seconds** as it can freeze too deeply and cause damage to the skin, resulting in scars and/or nerve damage. Freezing the skin tag too deeply can lead to local tissue death.

F. WARNINGS

KEEP OUT OF REACH OF CHILDREN.



Danger - Extremely Flammable Aerosol.

Pressurised container: May burst if heated. Keep away from heat, hot surfaces, sparks, open flames and other ignition sources. No smoking. Do not pierce or burn, even after use. Do not spray on open flame or other ignition source. Protect from sunlight. Do not expose to temperatures exceeding 50°C/122°F.

- **Avoid Contact with Eyes.** Contact with eyes may cause blindness. Do not rub or touch eyes if Cryotag has been applied to hands. If product gets into eyes, flush with water for 15 minutes and get medical help right away.
- **Avoid Inhaling and Use Only in Well-Ventilated Areas.** It is uncertain how the contents will affect the lungs. Breathing the contents is hazardous.
- **Do Not Eat.** Contents are toxic and may cause serious internal damage. If eaten by accident, get medical help or contact a poison control center immediately.
- **For External Use Only.**
- **Use ONLY as Directed.** Intentional misuse by deliberately concentrating and inhaling contents can be harmful or fatal. Misuse of Cryotag may result in burns and permanent scarring of healthy tissue or blindness. Failure to follow instructions may result in serious injury or harm.



CryoConcepts, LP
205 Webster Street,
Bethlehem, PA 18015 USA

EC REP
Emergo Europe
Prinsessegracht 20,
2514 AP The Hague,
The Netherlands



Distributed by Appia
Healthcare limited,
PO Box 1456, Guildford,
GU1 9RU, United Kingdom.

G. POSSIBLE RISKS AND SIDE EFFECTS

Freezing destroys the tissue and begins the removal process. As a result of freezing, a blister may form (sometimes containing blood) within a few days.

- A stinging, itching or aching feeling may occur during or after the freeze treatment.
- Freezing treatments may cause a change in skin color. Dark skin may become lighter in appearance or light skin may become red or darker in appearance.

H. HOW TO USE CRYOTAG?

Cryotag is very easy to use. Please follow the 4 simple steps below.

- **Read all instructions carefully before using. DO NOT SPRAY DIRECTLY ONTO SKIN. DO NOT APPLY THE APPLICATOR IMMEDIATELY. Do not hold the spray can near your face or over parts of your body or clothing. Do not let any liquid make contact with your skin.**

Getting ready: you will need a clock or household timer that counts time in seconds as a way to time your treatment.



Hold a Cryotag precision-tip applicator by the long plastic tube/stick. Attach the precision-tip applicator stick to the cap by firmly twisting it into the opening on bottom of the cap as shown.



Place the canister on a table or other sturdy surface. This makes it easier to press the cap. Insert the precision-tip applicator with the cap into the shield on the canister. Keeping the canister upright, press the cap all the way down for 3 seconds. You will hear a hissing sound.

For any other question or information, please visit www.cryotag.co.uk

©2019 CryoConcepts, LP Bethlehem, PA 18015
This product is covered by one or more patents.
See www.CryoConcepts.com/patents



Remove the precision-tip applicator and the cap from the shield and the canister. Hold with the tip facing down. Do not apply the applicator immediately; wait 15 seconds before applying to the skin tag. This permits the applicator tip to reach the correct temperature.



Grasp the skin tag with the tweezers and gently pull it away from the base. Place the tip of the applicator to the stalk of the skin tag that is to be frozen. Apply the precision-tip applicator to the stalk continuously for 40 seconds and no longer. Do not dab.

Note: It is normal to feel an aching, stinging sensation during treatment and for a short time afterwards. Remove the precision-tip applicator from the cap and discard the precision-tip applicator after every single usage. No special method of disposal is required. Retain the canister, cap, remaining precision tip applicators, instruction leaflet in original packaging for future use.

I. FOLLOW-UP TREATMENT

Skin tags may be removed in as little as 1 treatment. If the skin tag has not fallen off in 2 weeks, you may treat the skin tag again, up to 4 treatment cycles. If the skin tag persists, stop using the product and talk to your doctor.

- **Be alert to signs of a developing infection:**
 - Redness or hardening of surrounding area.
 - Treated area is hot to touch.
 - Yellow or green discharge from treated area.
 - Change in odor of discharge.
 - Increasing or unusual pain in treated area.
 - Fever.
 - Painful red streaks below skin surface.
- **To help you avoid infection:**
 - Keep the treated area clean.
 - Protect blisters, if necessary, with a bandage; do not puncture blisters.
 - Do not pick at, or scratch the treated area.
- **To help with any discomfort:** If your skin tag is rubbing against some clothes or jewellery after the treatment and you feel some discomfort, you can protect it by applying a plaster.
- Swimming and showering are permitted.
- Talk to your doctor if you have burning, aching or stinging sensations that do not go away or if you think the treated area is getting infected.